BRAUN


Sensian™ Swipe 5

Temple swipe thermometer

BST200

GB	English	2
DE	Deutsch	9
ES	Español.....	17
FR	Français	24
IT	Italiano	32
NL	Nederlands	39
PT	Português	46

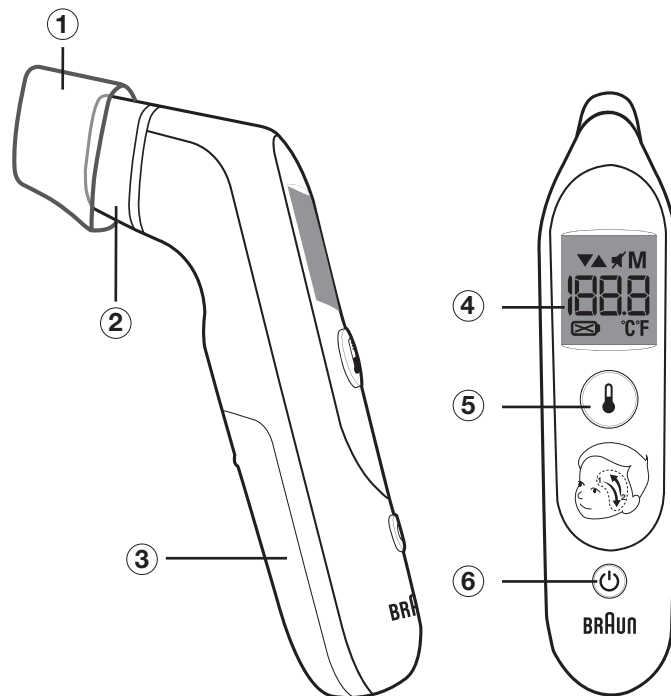 Kaz Europe Sàrl
A Helen of Troy Company
Q-Center, Chaux 4
CH-1030 Bussigny - Switzerland

© 2023, All rights reserved.

www.BraunHealthcare.com/uk_en/

Certain trademarks used under license from
The Procter & Gamble Company or its affiliates.

Sensian™ and TempleSwipe™ Technology are
trademarks of Helen of Troy Limited and/or its
affiliates.

www.BraunHealthcare.com/us_en/patents

Made and printed
in China

MD **REF** BST200WE

P/N: A004712R1
10JAN23

Product description (See page 1)

- ① Protective cap
- ② Infrared temperature scanner
- ③ Battery door
- ④ LCD screen
- ⑤ Temperature button
- ⑥ Power button

Thank you for purchasing the Sensian™ Swipe 5 Temple swipe thermometer (BST200). This thermometer is a high quality product incorporating the latest technology and tested in accordance with international standards. With its unique technology, the product can provide a stable, heat-interference-free reading with each measurement. The instrument performs a self-test every time it is switched on to help provide accurate measurements.

Please read these instructions carefully before using this product and keep the instructions and the thermometer in a safe place.

Use of this thermometer is not intended as a substitute for consultation with your doctor.

Intended use and contraindications

This Braun thermometer is a handheld, battery-powered, infrared thermometer intended for the intermittent measurement of human body temperature in a home-use environment on people of all ages, except pre-term and small-for-gestational-age babies.

These thermometers are not intended for clinical use in a professional environment and are for home usage only. This thermometer is not intended to interpret hypothermic temperatures.

Warnings and precautions

Keep out of reach of children under 12 years.

NEVER use the thermometer for purposes other than those it has been intended for. Please follow the general safety precautions when using on children.

NEVER immerse the thermometer into water or other liquids (not waterproof). For cleaning and disinfecting please follow the instructions in the "Care and cleaning" section.

DO NOT store this thermometer in temperature extremes below -25 °C or over 55 °C (below -13 °F or over 131 °F) or in excessive humidity (above 95% non-condensing relative humidity).

If thermometer is stored in a location that is cooler or warmer than where it is being used, let it sit in the patient's room for 10 minutes before taking the measurement.

DO NOT use the thermometer if there are signs of damage on the scanner or on the thermometer itself. If damaged, **DO NOT** attempt to repair the product.

NEVER scratch the sensor with a sharp object when cleaning.

This thermometer consists of high-quality precision parts. **DO NOT** drop the instrument. Protect it from severe impact and shock. **DO NOT** twist the instrument or the measuring sensor.

Temperature elevation may signal a serious illness, especially in neonates and infants, or in adults who are old, frail, or have a weakened immune system. Please seek professional advice immediately when there is a temperature elevation and if you are taking temperature on:

- Neonates and infants under 3 months (Consult your doctor immediately if the temperature exceeds 37.4 °C [99.4 °F])
- Patients over 60 years of age (Fever may be blunted or absent in older patients)
- Patients having diabetes mellitus or a weakened immune system (e.g., HIV positive, cancer chemotherapy, chronic steroid treatment, splenectomy)
- Patients who are bedridden (e.g., nursing home patient, stroke, chronic illness, recovering from surgery)
- A transplant patient (e.g., liver, heart, lung, kidney)

DO NOT allow children to take their temperatures unattended.

Please consult your doctor if you see symptoms such as unexplained irritability, vomiting, diarrhea,

dehydration, changes in appetite or activity, seizure, muscle pain, shivering, stiff neck, pain when urinating, etc., even in the absence of fever.

Even in the absence of fever, those who exhibit a normal temperature may still need to receive medical attention. People who are on antibiotics, analgesics, or antipyretics should not be assessed solely on temperature readings to determine the severity of their illness.

DO NOT modify this equipment without the authorization of the manufacturer.

Why Braun Sensian™ Swipe 5?

Measurement in 3 seconds

The Braun TempleSwipe™ Technology allows measurement of temple temperature in just 3 seconds.

Accurate and reliable

Thanks to Braun TempleSwipe™ Technology, the Sensian™ Swipe 5 thermometer captures the heat naturally given off by the temple to calculate oral equivalent temperature.

Easy to use

The Sensian™ Swipe 5 thermometer is non-invasive. A measurement can be taken even while a child is sleeping.

Design is easy to hold and use.

The Braun Sensian™ Swipe 5 thermometer is safer to use on a child than a rectal thermometer and is easier to use than other methods.

Safe and hygienic

Helps minimize spreading of germs.

Completely safe for use on children and adults.

Auto-Display memory

The last reading is automatically displayed for 2 seconds when the unit is switched ON.

How does Braun Sensian™ Swipe 5 work?

The Braun Sensian™ Swipe 5 thermometer measures infrared energy radiated from the skin above the eyebrow area, the temple, and the surrounding tissue. Temperature readings obtained by scanning the temple area will provide the greatest accuracy.

This captured energy is collected through the sensor and converted to an oral equivalent temperature value.

The Braun Sensian™ Swipe 5 thermometer has been clinically tested and proven to be safe and accurate when used in accordance with its operating instruction manual.

Normal temperature is a range. The range of normal also varies from person to person and fluctuates throughout the day. It is, therefore, important to determine your normal temperature range.

Colour-coded fever guidance


Colour-coded fever guidance helps you to better understand the meaning of your child's temperature with the color indicated on the display. The screen displays green for a normal temperature, yellow for a fever, and red for a high fever.

Color range	Reading	Meaning
Green	≥ 35.8 – ≤ 37.4 °C (≥ 96.4 – ≤ 99.4 °F)	Normal
Yellow	> 37.4 – ≤ 38.5 °C (> 99.4 – ≤ 101.3 °F)	Fever
Red	> 38.5 °C (> 101.3 °F)	High fever

How to use your Braun Sensian™ Swipe 5

1. Power on

Remove protective cap. Press and release the power button once. Backlight will come on and the warm up sequence starts. The last reading is automatically displayed for 2 seconds.

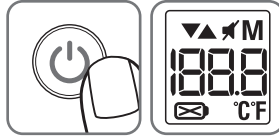When device displays STEADY "00.0" on the screen, it is ready to take a measurement.

2. Position

Place the thermometer against the skin, on the forehead just above the eyebrow. Be sure to hold the thermometer flush against the forehead.

NOTE: ALWAYS remove the protective cap before taking a temperature reading.

3. Take temperature


Press the temperature button and slowly swipe down forehead to the ear and **BACK AGAIN** while the short beeps are heard. Be sure to hold the thermometer flush against the forehead for the duration of the measurement.

When the temperature button is pressed a dashed line will appear on the screen and a soft beeping is heard to let you know the thermometer is working. **DO NOT** remove the thermometer from the skin until you hear the long, 2-second confirmation beep.

4. Read temperature

At the confirmation beep, remove the thermometer from the skin and read the temperature. The screen will illuminate the appropriate color for the temperature reading. Green for a normal temperature, yellow for a fever, and red for a high fever.

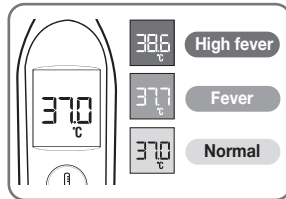5. To turn off

Press the power button. To conserve battery power, the unit will automatically shut down after 60 seconds of no use.

Temperature taking hints

It is important to know each individual's normal temperature when they are well. This is the only way to accurately diagnose a fever. Take multiple readings when healthy to determine normal temperature.

A child's normal temperature can be as high as 37.7 °C (99.9 °F) or as low as 36.1 °C (97.0 °F). Be sure to note this unit reads 0.5 °C (0.9 °F) lower than a rectal digital measurement.

An individual must be inside for 30 minutes before taking a measurement.

Hold the thermometer flush against the forehead and swipe slowly from above the eyebrow down to the temple, and back again. Swiping too slow or too fast may impact your reading.

Swipe down the temple to the ear and back again. Use the beeps to help pace your swiping so that you are in the starting position when the long beep is heard.

DO NOT remove the thermometer from the forehead before hearing the final beep.

Nursing mothers should not take a temperature measurement on themselves or have their temperature measured by another individual while or immediately after nursing a baby, as nursing can affect body temperature.

DO NOT use the Sensian™ Swipe 5 thermometer in high humidity environments.

Individuals should not drink, eat, or be physically active before/while taking the measurement.

Before placing the thermometer sensor onto the forehead, remove dirt, hair, or sweat.

Use an alcohol swab to carefully clean the sensor and wait for 5 minutes before taking a measurement on another Individual. Wiping the forehead with a warm or cool cloth may impact your reading. It is advised to wait 15 minutes before taking a reading.

ALWAYS take the temperature in the same location, since the temperature readings may vary from different locations on the forehead.

In the following situations, it is recommended that three temperatures be taken and the highest one used as the reading:

- Infants within the first 100 days of their birth.
- Children under three years of age with a compromised immune system and for whom the presence or absence of fever is critical.
- When the user is learning how to use the thermometer for the first time until he/she has familiarized himself/herself with the instrument and obtains consistent readings.


Last reading call

The Braun Sensian™ Swipe 5 thermometer includes a memory feature to recall the last temperature reading. The last temperature taken is stored in its memory and will be automatically displayed for 2 seconds when it is turned on.

How to change temperature between °C and °F and changing the sound mode

This thermometer has the ability to switch temperature units from Celsius to Fahrenheit as well the ability to mute the device. When the device is in mute mode, all audible sounds will be disabled. The mute icon will be displayed on the LCD at all times when the device is in mute mode.



1. While the device is off, press and hold the  button for 6 seconds until SET appears on the screen.
2. Release the measurement button  when the desired temperature scale and mute option are shown. The selected icons will flash on the screen to confirm the new setting, then the thermometer will turn off automatically.

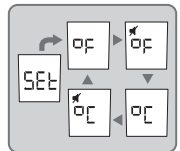Calibration

This product has been calibrated at the time of manufacture. If used according to the use instructions, periodic calibration is not required. This device is not intended as a substitute for regular check-ups by your doctor, please consult your doctor if you have any doubt about the temperature reading.

Errors and troubleshooting

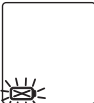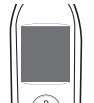Error message	Situation	Solution
	Error symbol: during taking measurement A "Err" sign appears when the temperature measured is BELOW 34 °C or 93.2 °F. The mark "Err" appears if the operating temperature is ABOVE 42.2 °C or 108 °F.	When you receive this symbol, press the power button to reset the thermometer before taking a new measurement.
	When ambient temperature is too high or too low to take an accurate reading, device will provide an "Err" message with flashing up/down arrow icon. The error message will appear after start up and memory recall.	When you receive this symbol, press the power button to reset the thermometer before taking a new measurement. If the ambient temperature is not within the allowed range of 15 °C - 40 °C (59 °F - 104.0 °F), you will not be able to take a reading.
		
	Low battery When 20% of battery life is left, the display FLASHES the low battery warning symbol, however, the device can still work until the battery life has 0% left.	Replace batteries. See instructions on page 7.
	Dead battery indication If the flashing battery icon is the only symbol shown on the display, the device cannot work.	Replace batteries. See instructions on page 7.
	Blank display Thermometer does not have power.	Please check if the batteries have been loaded correctly. Also check polarity (<+> and <->) of batteries.
	All segments display or Blank display System error.	Replace batteries. See instructions on page 7. If error persists, call customer service.


Care and cleaning

Use an alcohol swab or cotton swab moistened with alcohol (70% Isopropyl) to clean the thermometer casing and the measuring probe. Ensure that no liquid enters the interior of the thermometer. **NEVER** use abrasive cleaning agents, thinners or benzene for cleaning and **NEVER** immerse the instrument in water or other cleaning liquids. Take care not to scratch the surface of the lens by replacing the protective cap after use.

Replacing the batteries

The Sensian™ Swipe 5 thermometer comes with 2 AAA batteries.

To install the batteries, slide open the battery door as shown. Insert 2 AAA batteries, correctly positioning with positive "+" and negative "-" as shown. Replace the battery door and close securely.

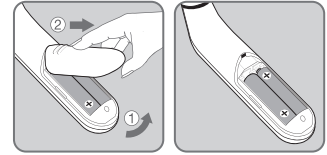Replace the batteries when the flashing battery symbol appears on the LCD screen. Dispose of batteries according to local waste management directions.

Only discard empty batteries. They should not be disposed of in household waste, but at appropriate collection sites or at your retailer.

2 Year Limited guarantee

Please read all instructions before attempting to use this device. Please retain the receipt as proof of and date of purchase. The receipt must be presented when making any claim within the relevant guarantee period. Any claim under the guarantee will not be valid without a proof of purchase. Your device is guaranteed for Two Years (2 Years) from date of purchase. This guarantee covers defects in materials or workmanship that occur under normal use; defective devices meeting these criteria will be replaced free of charge.

The guarantee DOES NOT cover defects or damage resulting from abuse or failure to follow the user instructions. The guarantee becomes void if the device is opened, tampered with, or used with non-Braun branded parts or accessories, or if repairs are undertaken by unauthorized persons.

Accessories and consumables are excluded from any warranty.

For support requests, please visit www.BraunHealthcare.com/uk_en/ or find contact information at the end of this owner's manual.

This Guarantee is applicable to Europe, Russia, Middle East and Africa only.

UK Only: This does not affect your consumer statutory rights.

Inquiry

In Europe, any serious incident (e.g. death, life-threatening injury, surgical intervention, ...) that has occurred in relation to this thermometer should be reported to Kaz Europe (see Helen of Troy Consumer Service Numbers) and to the competent authority of your member state.

REF, LOT and Serial numbers of your product must be provided as they are essential to record and follow up your inquiry or claim. The LOT and Serial numbers are located in the battery compartment. The production date is given by the LOT number and can be deciphered as explained below:

The first 3 numeric digits in the LOT number represent the day of the year of manufacture. The next 2 numeric digits represent the last two numbers of the calendar year of manufacture and the letter(s) at the end designates the production site. (E.g.: LOT No.: 12313tav this product was made on the day 123, year 2013 at production site "tav")

Product specifications


Type: Temporal swipe (BST200)
Measuring range: 34 °C – 42.2 °C (93.2 °F – 108 °F)
Resolution: 0.1 °C (0.1 °F)
Accuracy (Laboratory) ±0.2 °C for 35.0 °C to 42.0 °C
(±0.4 °F for 95.0 °F to 107.6 °F)

Clinical accuracy and procedures ±0.3 °C (±0.5 °F) outside of this range are available upon request.

Display: Liquid Crystal Display, 4 digits plus special icons
Acoustic: Audio:

Single long beep for 2 sec duration
Operating temperature: 15 °C – 40 °C (59 °F – 104.0 °F)

Automatic Switch-off: Approx. 60 seconds after last measurement has been taken
 Weight: 92g (±5g) (with batteries), 69g (±5g) (w/o batteries)
 Long term storage ranges
 Storage/transport temperature: -25 °C to 55 °C (-13 °F to 131 °F)
 Humidity: 15–95% non condensing
 Battery: (2) AAA Batteries - at least 1000 measurements
 Pressure: 700–1060 hPa (0.7–1.06 atm)
 This thermometer is an adjusted mode thermometer that provides an oral equivalent reading.

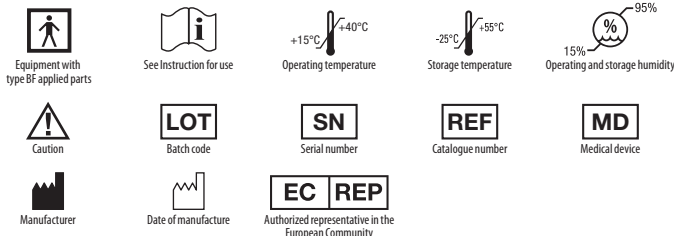Internally Powered Equipment

Continuous Operation

IP22: Protected against solid foreign objects of 12.5 mm diameter and greater. Protected against vertically falling water drops when the device is tilted up to 15 degrees.

MEDICAL ELECTRICAL EQUIPMENT needs special precautions regarding EMC. NOTE: Do not use this device in the presence of electromagnetic or other interference outside the normal range specified in EN 60601-1-2. For detailed description of EMC requirements, please visit website (Support / Document center) or contact service center. Portable and mobile RF communications equipment can affect MEDICAL ELECTRICAL EQUIPMENT.

Clinical accuracy and procedures are available upon request.

This medical device carries the CE mark and is manufactured in conformity with RoHS Directive 2011/65/EU and other applicable directive and/or regulations as specified in the EU Declaration of Conformity.

E 1965-98 for the thermometer system. Full responsibility for the conformance of the product to the standard is assumed by Kaz Europe Sarl, Switzerland.

ASTM laboratory accuracy requirements in the display range of 37 to 39°C (98 to 102°F) for IR thermometers is ±0.2°C (±0.4°F), whereas for mercury-in-glass and electronic thermometers, the requirement per ASTM Standards E667-86 and E1112-86 is ±0.1°C (±0.2°F).

 Please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at your local retailer or at appropriate collection points provided in your country.

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- ① Schutzkappe
- ② Infrarot-Temperatur-Scanner
- ③ Batteriefach-Deckel
- ④ LCD-Bildschirm
- ⑤ Temperaturtaste
- ⑥ Ein/Aus-Taste

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Helen of Troy Consumer Service Numbers

Algeria

☎ +213 34 18 71 42

Austria

☎ +43 (0)1 360 277 1225

Bahrain

☎ +973 17228887

Belgium/Luxemburg

☎ +32 (0) 2 620 01 01

Bosnia

☎ +387 (33) 902 911

Bulgaria

☎ +359 2 492 87 05

Croatia

☎ +385 (01) 3444 856

Cyprus

☎ +35722337593

Cyprus Northern region

☎ +90 (0) 392 22 91 461

Czech Republic

☎ +420 (0) 2 25 43 97 69

Denmark

☎ +45 35 15 80 40

Finland

☎ +358 (0)9 81 71 00 14

France

☎ +33 (0) 1 85 14 80 95

Germany

☎ +49 (0) 21 173 749 003

Greece

☎ +30 211 180 94 56

Hungary

☎ +36 (06) 1 429 2216

Iceland

☎ +354553100

Ireland (The Republic Of)

☎ +353 (0) 1 525 1809

Israel

☎ +1 800250221

Italy

☎ +39 02 3859 1183

Jordan

☎ +962 64168003

Kuwait

☎ +965 184 0011

Lebanon

☎ +961 (01) 512002

Netherlands

☎ +31 (0) 78 201 8001

Norway

☎ +47 23 50 01 20

Poland

☎ +48 22 512 39 02

Portugal

☎ +351 21 060 8045

Qatar

☎ +974 4 4075048

Romania

☎ +40 213 520 167

Saudi-Arabia

☎ +966 920001414

Slovakia

☎ +421 2 5011 2115

Slovenia

☎ +386 (01) 888 86 74

South Africa

☎ ++27 (0) 11 844 6190

Spain

☎ +34 913 754 176

Sweden

☎ +46 (0)8 5199 3097

Switzerland

☎ +41 (0) 22 567 5200

Tunisia

☎ +216 98 780237

Turkey

☎ +90 (0) 212 217 66 60

UAE

☎ +971 4 349 4014

United Kingdom/ Northern Ireland

☎ +44 (0) 207 949 0115

**Helen
of Troy**

**Creative Department
Artwork Specifications**

Brand: Braun
Category: Thermometer
Model: **BST200WE**
Artwork Part #: A004712R1
Die Line Part #: NA
Subject: Owners Manual
Region: EMEA
Flat Size: **MM: W215.9 x H177.8**
IN: W8.5 x H7.0
Folded Size: **MM: W107.95 x H177.8**
IN: W4.25 x H7.0
Scale: 1/1
Material: 120g wood free paper
(uncoated) for Cover
80g wood free paper
(uncoated) for Interior
Page count: 54
Revision: 3
Date: **10JAN23**
Release Date: 13JAN23
Rerelease Date:

Colors:

Dielines (Do not print)

Spot Colors

Cyan
0%

Magenta
0%

Yellow
0%

Black
100%

Special Instructions:

Quality Requirement of Artwork and Quality
Clarification Process of Artwork Printing:
Meet Eng-QS-06&02

Helen of Troy
Creative Services
Marlborough, MA 01752 USA
+1 508 490 7000